


PTR-QMS SERIES



THE WORLD'S
LEADING
PTR-MS
COMPANY PRESENTS

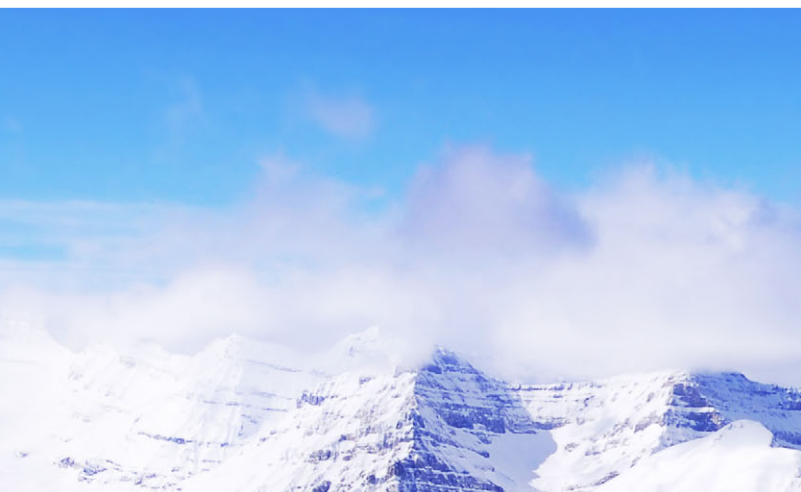

PTR-QMS 500

LoD < 30 pptv

The IONICON PTR-QMS 500 instrument is a **very sensitive detector** for volatile organic compounds (VOCs) that allows for **continuous VOC quantification**.

Our PTR-QMS series product, based on quadrupole mass spectrometry, combines **low online detection limits** reaching two-digit pptv-levels, with high selectivity and a very fast response time.

Direct injection of sample gases **without preparation** contributes to the **speed and simplicity** that is common to all our instruments.

Selective qualitative and quantitative analysis of trace compounds in remarkably short measurement times combined with a system **linearity range** covering over **six orders of magnitude** satisfy the highest claims of the world's best scientists.

Our unique **soft ionization (PTR) technology** together with our extensive experience in gas-phase ion chemistry and engineering of scientific instruments accounts for the **reliability, very low detection limit, very low mass fragmentation, fast response time and robustness** of our PTR-MS systems.

Mass Range 1-512 amu

- > Low detection limit < 30 pptv
- > Real-time VOC quantification
- > Soft & efficient ionization technology
- > Fast and easy direct sample injection

Find out more:

www.ionicon.com/products

PTR-QMS 500



IONICON PTR-QMS 500 SPECIFICATIONS*

- Mass range: 1-512 amu
- Resolution**: < 1 amu
- Response time: 100 ms
- Sensitivity (Benzene): > 80 cps/ppbv
- Detection limit**: < 30 pptv
- Linearity range**: 30 pptv - 10 ppmv
- Adjustable flow: 80 - 800 sccm (at 60°C)
- Inlet system (Different inlet systems available on request):
 - 1.2 m long inlet hose - with internal inert (PEEK) capillary
 - Inlet system heating: up to 180°C (356°F)
- Reaction chamber heating range: 45 - 120°C (104 - 248°F)
- Power supply and max. consumption: 115/230 V, 1500 W
- Dimensions (w x h x d): 60x91x60 cm (23.7x35.9x23.7 in.)
- Weight (incl. SRI): 129 kg (285 lbs)
- Main interfaces:
 - 1x Touch screen display
 - 1x Ethernet 10/100 Mbit RJ45 (TCP/IP)
 - 8x DO, 4x DI, 4x AI, 4x AO

*Specifications are subject to change without prior notice.
 Product pictures and illustrations may differ from actual configuration.
 **Detection limit, linearity range and resolution are dependent on the substances measured, integration time and system set-up.

TECHNOLOGY

The innovative technology IONICON PTR-MS products are based on is Proton Transfer Reaction - Mass Spectrometry (PTR-MS).

This unique soft ionization is realized by proton transfer from H_3O^+ ions to all compounds with a higher proton affinity (PA) than water. Common constituents of air such as N_2 , O_2 , Ar, CO_2 etc. have lower PAs than H_2O and are therefore not ionized. This is one of the main reasons for our market-leading low online detection limit for trace compounds and due to precisely controlled ion source and drift tube parameters, absolute quantification of VOC concentrations is possible.

PTR+SRI-MS

SWITCHABLE REAGENT IONS

IONICON invented the Switchable Reagent Ions (SRI) technology featuring NO^+ and O_2^+ as additional precursor ions created in the new IONICON ULTRA-PURE ion source. The benefits are extraordinary as not only isomeric VOCs can be separated and instantaneously quantified but also substances with a smaller PA than the PA of H_2O can now be detected with PTR-MS instruments.

4th GENERATION PTR-MS

The PTR-QMS 500 is incarnating the technological evolution of more than a decades's experience in PTR-MS we have at IONICON.

With the 4th generation of IONICON PTR-MS instruments, featuring technology such as the new ULTRA-PURE ion source and many other refinements, the new PTR-QMS 500 is outpacing most offerings in the market.

The PTR-QMS 500 is the rational choice for users in the need of a very sensitive real-time mass spectrometer for VOC analysis. It has all the benefits of the new PTR-QMS 500 series but a lower price tag than the High-Sensitivity PTR-QMS 500 instrument, the series flagship model.

The PTR-QMS 500 represents the 4th and most advanced generation of IONICON PTR-MS systems. It is the successor of the IONICON Standard PTR-MS, of which various versions and generations are in use all over the world.

The PTR-QMS 500 can be upgraded to a High-Sensitivity PTR-QMS 500 system if the noteworthy detection limit of < 30 pptv is not sufficient but the market-leading < 1 pptv of the High-Sensitivity PTR-QMS 500 is needed.

We provide flexible inlet solutions such as the integrated multi-position capability, calibration solutions and the possibility to optionally upgrade the PTR-QMS 500 series with multiplexing and data interface units thus bringing down the costs per sampling point and simplifying the integration in existing monitoring network infrastructure.